Smart**W**atcher®
SPIRIT

USER GUIDE
BENUTZERHANDBUCH
MODE D'EMPLOI
GUIDA UTENTE
GEBRUIKERSHANDLEIDING

ABOUT THIS WATCH

The Smartwatcher emergency call watch allows you to speak with selected friends and family or with care professionals in order to arrange assistance.

Optimised for safety, this watch only works together with the Smartwatcher subscription service of your choice. Other features are not available on this device. In daily use, the watch shows the time.

Find more information on
www.smartwatcher.com



**Welcome to the
Smartwatcher experience!**

CONTENTS

1. What's in the box.....	5
2. Watch set up.....	6
Step 1 - Charge your watch	
Step 2 - Register your watch	
3. How to use your watch	7
• Power on	
• Power off	
• Restart	
• Charging	
• To see the time	
• To trigger an alert	
• To end an alert	
• The Watch Management Portal	
• Test your watch	
4. How to take care of your watch.....	15

1. WHAT'S IN THE BOX

- Watch
 - Watch face (1)
 - Alarm button (2)
 - Power on/off button (3)
 - Microphone (4)
 - Speaker (5)
- Charging cable & cradle (6)
- USB power adapter (7)
- Quick starter,
user guide,
warranty

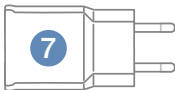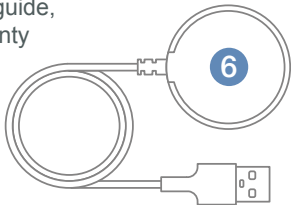2. WATCH SET UP

Step 1 - Charge your watch

Fully charge your watch before using it for the first time.

1. Plug the charging cable into the USB power adapter and connect it to the power socket.
2. Place the watch on the charging cradle. It takes about 2 hours to fully charge the watch.

Step 2 - Register and activate your watch online

1. Power on your watch by pressing and holding down the power on/off button **3** for 2-3 seconds.
2. Go to my.smartwatcher.com in your computer Internet browser.
3. Insert the watch's IMEI and ICCID numbers which are shown on the watch face **1**.
4. You can now follow the registration and activation process step by step.

TIP

When you add a helper's mobile phone number, that person will receive an SMS with a link to download and install our free Smartwatcher app. After installing the app, Helpers need to accept the invitation. Now, during an active alarm, they will be able to see the location of the person wearing the watch.

In order to manage your watch settings and the safety features please visit the Watch Management Portal at my.smartwatcher.com

3. HOW TO USE YOUR WATCH

Power on

Press and hold down the power on/off button **3** for 2-3 seconds.

Power off

1. Press the lower button **3** twice quickly.
2. When you see the „Power off“ watch face, please select:
 - a. the green icon (right)  to confirm you want to Power off the watch
 - b. the red icon (left)  to cancel the Power off

Restart

Press and hold down the power on/off button **3** for 8-9 seconds.

Charging

For charging instructions please refer to Step 1 of „Watch Set Up“.

To see the time

The watch face is activated by gently turning your wrist towards you. Alternatively, lightly press the power button **3** for 1 second.

To trigger an alert

Firmly press the large upper button **2** for 1-2 seconds until the watch vibrates and you hear a beep tone. The watch face **1** now shows a 5 second countdown, which allows you to end a false alarm. If you do not end the alarm, a red alarm symbol appears on the watch face **1** and the alarm goes out to your helpers.

Useful information

- Please plan to charge your watch at least once each day.
- Use only the originally-supplied charging set. Several factors can impact how long the battery will last between charging:
 - The strength of your mobile network connection
 - Intensity of use
 - Extreme temperatures (very warm or cold)

IMPORTANT!

Network coverage is not available everywhere. You will not be able to trigger an alarm when cellular service is not available. If there is no mobile network connection, the watch will show the following indicator: .

To end an alert

An active alarm can be ended in three different ways:

1. Automatically, when the call ends. The alarm will automatically end after the last helper hangs up the phone (if the call lasted at least 10 seconds).
2. Manually by the person wearing the watch. You can end the alarm on your watch by tapping the red alarm icon on the watch face. Then tap the green confirmation symbol on the right side. This returns you to the watch face showing the time.
3. Online, in the Watch Management Portal. The watch administrator can end an alarm by navigating into the „Alarm History“ menu and clicking on „End alarm“.

The Watch Management Portal

In the Watch Management Portal (my.smartwatcher.com) the watch administrator can manage watch settings, helpers, personal information and safety features.

Test your Watch

EN

IMPORTANT!

To support the performance of your watch, please **test the connectivity at home and any areas where you will regularly wear the watch.** While the SPIRIT/SPIRIT Petite shows excellent cellular network connectivity, low network connectivity in general can impact the watch's ability to send an alarm.

Please note that if there is no cellular connectivity on the watch, a call cannot be established.

NOTE

With SPIRIT/SPIRIT Petite specifically, the watch works very hard to connect. When cellular connectivity in the area is low, it may become warm and give itself an automatic cooling period. When cooling, the watch shows that it is in airplane mode and cannot connect a call. After cooling, the watch will return to the watch face visual or give you the option to press „reload.“

We recommend that you test your watch alarm periodically (i.e. once a month) to ensure it is connecting properly.

4. HOW TO TAKE CARE OF YOUR WATCH

- Always clean your Smartwatcher SPIRIT/SPIRIT Petite with a clean, soft cloth.
- If the device was exposed to clear water (fresh water), dry it carefully with a clean, soft cloth. If the device was exposed to other liquids, such as salt water, swimming pool water, soapy water, oil, perfume, sunscreen, hand cleanser, or chemical products such as cosmetics, contrary to an existing IP classification, which only applies to clear fresh water, you should rinse it with fresh water and dry it carefully with a clean, soft cloth. Contact with other liquids can affect both the performance and appearance of the device.

For further instructions or help please visit

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© 2021 Smartwatcher Technologies AG
Smartwatcher Technologies AG